

ARMY RESEARCH LABORATORY



The Detection of Doped Cr Presence in $\text{Sr}_{0.6}\text{Ba}_{0.4}\text{Nb}_2\text{O}_3$ by SIMS

Unchul Lee

ARL-TN-132

May 1999

19990526 094

Approved for public release; distribution unlimited.

DTIC QUALITY INSPECTED 1

Abstract

Photorefractive materials like SrBaNb (SBN) have optical properties that can have several military applications. To improve the figures of merit for these materials, it is very important to have better control of Cr dopant in the SBN. Since the amount of Cr dopant is so small, the usual analytical techniques are not adequate to analyze this material. In this report, the incorporation of Cr into SBN using secondary ion mass spectrometry (SIMS) measurement with an electron gun is confirmed.

Contents

1. Introduction	1
2. Experiment	2
2.1 <i>Equipment</i>	2
2.2 <i>Depth Profiling</i>	3
3. Results and Discussion	4
3.1 <i>Presence of Cr</i>	4
3.2 <i>Uniformity</i>	4
4. Summary	6
References	6
Distribution	7
Report Documentation Page	9

Figures

1. Secondary ion mass spectrometry (SIMS)	2
2. SIMS data for nominal value samples of 0.005%, 0.010%, and 0.025% Cr	5

Tables

1. Average and normalized counts for three nominally doped Cr samples	4
2. Counts for spots A, B, and C for three nominally doped Cr samples	5

1. Introduction

Photorefractive (PR) materials are used for optical beam fanning, beam coupling, optical storage, and optical image processing [1]. These optical applications are critical for military purposes, such as protecting optical sensors. To maximize the protection capability, it is very important to increase the coupling coefficient and maintain a fast response time.

The control of the protection capability, i.e., the speed of response and coupling efficiency, requires a specific donor/acceptor ratio that is very difficult to realize using only a single dopant species. For this reason, the Rockwell Science Center introduced the double-doping concept and tried Cr^{3+} , which showed considerable promise in maintaining a fast response and a sufficiently high coupling coefficient.

SBN (Sr, Ba, Nb) crystals are grown along the (001) direction at the Rockwell Science Center using the Czochralski pulling technique from platinum crucibles in an oxygen atmosphere. To confirm the soundness of using Cr as a dopant, the Rockwell Science Center grew samples with different amounts of Cr dopant in the melts. In this study, we chose the materials doped with chromium only in the amounts of 0.005, 0.010, and 0.025 wt.% in the melts. The size of the samples was $6 \times 6 \times 6$ mm [2].

Analyzing these materials poses a challenge because the dopant concentration is too low for XPS/AES (x-ray photoelectron spectroscopy)/(auger electron spectroscopy), EDX/WDX (energy dispersive x-ray)/(wavelength dispersive x-ray) analytical techniques and because sample preparation for the AAS (atomic absorption spectroscopy) technique is very difficult. Therefore, we used the SIMS (secondary ion mass spectrometry) technique to prove the presence of Cr as a dopant as well as the amount by weight percent of Cr dopants. SIMS is the one of the most powerful techniques for the surface and depth profile analysis of materials.

SIMS typically provides detection limits between 1 ppm (part per million) and 1 ppb (part per billion), depending on instrumental transmission efficiency, detector efficiency, mass resolution, primary beam density, and matrix-related effects. SIMS instrumentation also provides a dynamic range of about nine orders of magnitude for secondary ion signal intensities.

This report details the dopant analyses conducted at the Army Research Laboratory of these photorefractive samples grown by the Rockwell Science Center.


2. Experiment

2.1 Equipment

The SIMS technique is widely used to investigate trace-level impurities and dopant distributions throughout the sample and to provide interface analysis for device fabrication. Many different kinds of SIMS exist, depending on the primary source and mass spectrometer, and pros and cons exist for each different type of SIMS instrument.

In this analysis, we used an ATOMIKA Ionprobe A-DIDA 3000 SIMS instrument with a quadrupole mass filter as a mass analyzer for the secondary ions (see fig. 1). The primary source is the magnetically confined hot cathode oscillating electron bombardment type. We extracted positive argon gas from the discharge using a two-level electrostatic immersion system and accelerated the primary ion beam to energies between 1 and 15 keV. In addition to this, the machine is equipped with cold cathode mode, which can produce O_2^+ ions as well. The source lens arrangement was designed as a telefocus system, providing a beam spot with current densities continuously variable between 10^{-10} and 10^{-4} A/cm² and with good homogeneity across a 5-mm-diam sample spot at a distance of about 55 cm from the source exit hole. We used a unipotential lens to focus the space-charge-broadened telefocus beam into a 100- μ m spot at the target surface. Two pairs of perpendicular deflection plates on the downstream side of the unipotential lens were used for beam scanning. Raster scanning of the primary beam, in conjunction with an electronic aperture, was used to reduce edge effects in depth profiling [3].

Figure 1. Secondary ion mass spectrometry (SIMS).

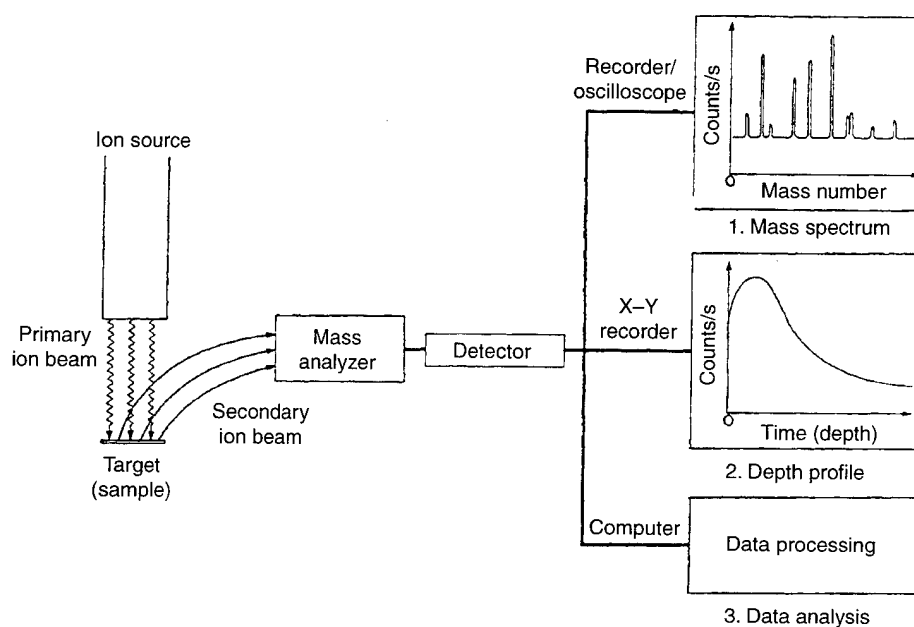

Proper selection of a narrow energy band pass for the secondary ions is essential for achieving good mass resolution in the quadrupole mass spectrometer. In the ATOMIKA SIMS instrument, this energy selection is performed in a modified parallel-plate condenser. Some of the secondary ions emitted from the sample enter the condenser through an entrance diaphragm mounted off-axis with respect to the condenser centerline and the quadrupole axis.

During SIMS analysis of samples with electrically insulating materials, a local charge builds up when samples are bombarded with a primary ion beam. This charge deflects the primary ion beam and accelerates the secondary ions in an undetermined manner. A simple but effective solution to this problem is to compensate this charge with an intense electron beam. Without charge compensation, no signal could have been obtained for the SBN samples analyzed in this study. Therefore, charge compensation was used for all the measurements we made during this study.

2.2 Depth Profiling

The three SBN samples were doped with Cr at nominal values of 0.005, 0.025, and 0.010 wt.%. For each sample at least three depth-profiling measurements were taken from different spots (A, B, and C). The primary oxygen ion beam energy was 6 keV with ~200 nA current on the sample surface. Since these samples were excellent insulators, it was not possible to collect any secondary ion signal because of the charging effect on the sample surface. When we bombarded the sample with a 1-keV electron beam of 2.7 A (filament), we then were able to obtain a measurable secondary ion signal. It would have been impossible to obtain the data without an electron gun.

3. Results and Discussion

3.1 Presence of Cr

With electron beam bombardment, we were able to obtain chromium and the major substrate constituents (Sr, Ba, Nb) ion signals from all three samples. For each depth profile, we took 50 layers of data points. Among the chromium isotopes, mass number 52 was chosen because Cr 52 is most abundant and there is no possibility of mass interference with other secondary ions.

We selected Sr counts as a benchmark. The average signal counts of Cr dopant were then normalized against Sr counts from the A spot of the 0.005 wt.% Cr-doped sample. As the data in table 1 show, these values are not proportional to the nominal values. One disadvantage of SIMS analysis is that it is not possible to get quantitative numbers unless standard samples with known dopant concentrations are available. It appears that the samples were either mislabeled or the gravimetric dopant quantity was incorrect.

3.2 Uniformity

To determine the uniformity of Cr dopant throughout the sample, at least three different spots (A, B, and C) were chosen from which to obtain depth-profile measurements. A difficulty was associated with the re-adjustment of the SIMS instrument whenever the measuring spot was changed. To make a reliable comparison, we needed a benchmark. For this reason, we included ions of three major constituents (Sr, Ba, Nb) as well as Cr counts for depth profiling. In each case, we normalized Cr counts against Sr counts for the A spot of the 0.005 wt.% sample, since the composition of the major constituents was identical for all the samples. Comparing the three graphs in figure 2, we note that the beginning of the 0.010 wt.% data points is different from the other two graphs. At the start of data collection for the SIMS measurement, there is an idiosyncrasy we must take into account, i.e., that the first few data points are always unstable. For this reason, the first few data points for figures 2(b) and (c) were not collected.

Average values were obtained from the compiled data. All the values were normalized against the 0.005 wt.% spot A Cr counts using Sr counts as a benchmark. From table 2, we conclude that there are variations in the uniformity of Cr dopant concentration throughout the surface and that the contents are not proportional as expected.

Table 1. Average and normalized counts for three nominally doped Cr samples.

Sample	Average counts	Normalized counts
0.005% Cr	50	50
0.010% Cr	70	32
0.025% Cr	150	114

Figure 2. SIMS data for nominal value samples of (a) 0.005%, (b) 0.010%, and (c) 0.025% Cr.

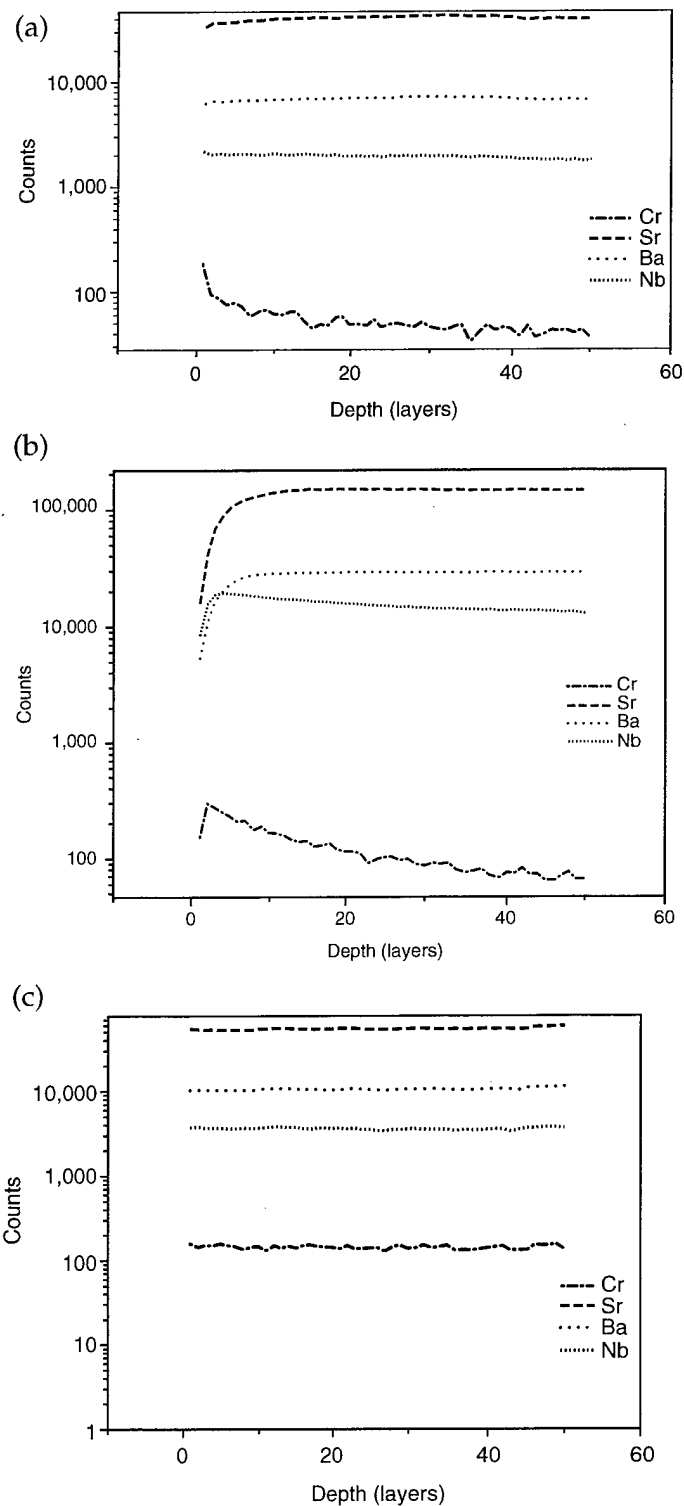


Table 2. Counts for spots A, B, and C for three nominally doped Cr samples.

Sample	Cr counts		
	Spot A	Spot B	Spot C
0.005%	50	30	30
0.010%	32	21	21
0.025%	114	105	210

4. Summary

First, we confirmed that Cr dopants had been incorporated into the SBN samples during the crystal growth process, a finding that was not possible using any other technique to date. We cannot say that the dopant proportionality is close to the proposed values. To calculate the absolute quantity of Cr dopant incorporated into the samples, we need standard samples that have exactly the same SBN matrix with a known quantity of dopants.

References

1. R. R. Neurgaonkar, W. K. Cory, J. R. Oliver, E. J. Sharp, M. J. Miller, W. W. Clark III, G. L. Wood, and G. J. Salamo, "Cr³⁺Sr_{0.6}Ba_{0.4}Nb₂O₆ Single Crystals for Photorefractive Applications," *Mat. Res. Bull.* **24**, 589 (1989).
2. R. R. Neurgaonkar, Technical Report SC5550.FR, Rockwell Science Center (January 1996).
3. A. Benninghoven, F. G. Rudenaur, and H. W. Werner, *Secondary Ion Mass Spectrometry*, John Wiley & Sons, New York (1987).

Distribution

Admnstr
Defns Techl Info Ctr
Attn DTIC-OCF
8725 John J Kingman Rd Ste 0944
FT Belvoir VA 22060-6218

Ofc of the Dir Rsrch and Engrg
Attn R Menz
Pentagon Rm 3E1089
Washington DC 20301-3080

Ofc of the Secy of Defns
Attn ODDRE (R&AT)
Attn ODDRE (R&AT) S Gontarek
The Pentagon
Washington DC 20301-3080

OSD
Attn OUSD(A&T)/ODDR&E(R) R J Trew
Washington DC 20301-7100

AMCOM MRDEC
Attn AMSMI-RD W C McCorkle
Redstone Arsenal AL 35898-5240

Army Rsrch Ofc
Attn AMXRO-GS Bach
PO Box 12211
Research Triangle Park NC 27709

CECOM
Attn PM GPS COL S Young
FT Monmouth NJ 07703

Dir for MANPRINT
Ofc of the Deputy Chief of Staff for Prsnl
Attn J Hiller
The Pentagon Rm 2C733
Washington DC 20301-0300

Hdqtrs Dept of the Army
Assist Secy of the Army for RD&A
Attn SARD-DOV H Fallin Rm 3E411
The Pentagon
Washington DC 20301-0103

Night Vsn & Elect Sensors Dirctr
Attn AMSEL-RD-NV-STD B Ahn
FT Belvoir VA 22060-5677

TECOM
Attn AMSTE-CL
Aberdeen Proving Ground MD 21005-5057

US Army Armament Rsrch Dev & Engrg Ctr
Attn AMSTA-AR-TD M Fisette
Bldg 1
Picatinny Arsenal NJ 07806-5000

Director
US Army CECOM Rsrch, Dev, & Engrg Ctr
FT Monmouth NJ 07703-5201

US Army Edgewood RDEC
Attn SCBRD-TD G Resnick
Aberdeen Proving Ground MD 21010-5423

US Army Info Sys Engrg Cmnd
Attn ASQB-OTD F Jenia
FT Huachuca AZ 85613-5300

US Army Natick RDEC Acting Techl Dir
Attn SSCNC-T P Brandler
Natick MA 01760-5002

Director
US Army Rsrch Ofc
4300 S Miami Blvd
Research Triangle Park NC 27709

US Army Simulation, Train, & Instrmntn
Cmnd
Attn J Stahl
12350 Research Parkway
Orlando FL 32826-3726

US Army Tank-Automtv Cmnd Rsrch, Dev, &
Engrg Ctr
Attn AMSTA-TA J Chapin
Warren MI 48397-5000

US Army Train & Doctrine Cmnd
Battle Lab Integration & Techl Dirctr
Attn ATCD-B J A Klevecz
FT Monroe VA 23651-5850

Distribution (cont'd)

US Military Academy
Mathematical Sci Ctr of Excellence
Attn MDN-A MAJ M D Phillips
Dept of Mathematical Sci Thayer Hall
West Point NY 10996-1786

Nav Surface Warfare Ctr
Attn Code B07 J Pennella
17320 Dahlgren Rd Bldg 1470 Rm 1101
Dahlgren VA 22448-5100

DARPA
Attn B Kaspar
3701 N Fairfax Dr
Arlington VA 22203-1714

Univ of West Virginia Physics Depart
Attn Abdul-Razzaq
Morgantown WV 26506

Hicks & Associates, Inc
Attn G Singley III
1710 Goodrich Dr Ste 1300
McLean VA 22102

Nat Inst of Stand & Techl
Attn P Amirtharaj
Semiconductor Bldg
Gaithersburg MD 20899

Palisades Inst for Rsrch Svc Inc
Attn E Carr
1745 Jefferson Davis Hwy Ste 500
Arlington VA 22202-3402

US Army Rsrch Lab
Attn AMSRL-DD J Rocchio
Attn AMSRL-CI-LL Techl Lib (3 copies)
Attn AMSRL-CS-AS Mail & Records Mgmt
Attn AMSRL-CS-EA-TP Techl Pub (3 copies)
Attn AMSRL-SE-EM U Lee (10 copies)
Adelphi MD 20783-1197

REPORT DOCUMENTATION PAGE			Form Approved OMB No. 0704-0188	
Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.				
1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE May 1999		3. REPORT TYPE AND DATES COVERED Final, 12/1/98 to 1/15/99
4. TITLE AND SUBTITLE The Detection of Doped Cr Presence in $\text{Sr}_{0.6}\text{Ba}_{0.4}\text{Nb}_2\text{O}_3$ by SIMS			5. FUNDING NUMBERS DA PR: AH94 PE: 62705A	
6. AUTHOR(S) Unchul Lee				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) U.S. Army Research Laboratory Attn: AMSRL-SE-EI email: ulee@arl.mil 2800 Powder Mill Road Adelphi, MD 20783-1197			8. PERFORMING ORGANIZATION REPORT NUMBER ARL-TN-132	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) U.S. Army Research Laboratory 2800 Powder Mill Road Adelphi, MD 20783-1197			10. SPONSORING/MONITORING AGENCY REPORT NUMBER	
11. SUPPLEMENTARY NOTES ARL PR: 9NE6B1 AMS code: 622705.H94				
12a. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release; distribution unlimited.			12b. DISTRIBUTION CODE	
13. ABSTRACT (Maximum 200 words) Photorefractive materials like SrBaNb (SBN) have optical properties that can have several military applications. To improve the figures of merit for these materials, it is very important to have better control of Cr dopant in the SBN. Since the amount of Cr dopant is so small, the usual analytical techniques are not adequate to analyze this material. In this report, the incorporation of Cr into SBN using secondary ion mass spectrometry (SIMS) measurement with an electron gun is confirmed.				
14. SUBJECT TERMS Cr dopant, SIMS, SBN			15. NUMBER OF PAGES 15	
			16. PRICE CODE	
17. SECURITY CLASSIFICATION OF REPORT Unclassified	18. SECURITY CLASSIFICATION OF THIS PAGE Unclassified	19. SECURITY CLASSIFICATION OF ABSTRACT Unclassified	20. LIMITATION OF ABSTRACT UL	